SyrVet™ Livestock Prod Instruction Manual

SyrVet, Inc. • Waukee, IA U.S.A. • Toll-Free: (800) 727-5203 • customerservice@syrvet.com • www.syrvet.com

⚠ WARNING

SAFETY INSTRUCTIONS:

THOROUGHLY READ AND FOLLOW ALL INSTRUCTIONS BEFORE USING THE SYRVET™ LIVESTOCK PROD.

- This prod is FOR USE ON APPROPRIATE LIVESTOCK ONLY.
- This prod is not intended to be used by children.
- Never use the prod on humans.
- KEEP PROD AWAY FROM CHILDREN.
- The prod produces a high voltage pulse. The electrodes should never be touched when the unit is on. Avoid touching the electrodes when batteries are installed in the unit.

- Clean and dry prod on a regular basis, removing batteries prior to cleaning and drying. Store prod in a dry location.
- KEEP PROD DRY. If prod gets wet, DISCONTINUE USE IMMEDIATELY and dry out to avoid corrosion and/or feedback. Do not use until prod is COMPLETELY dry.
- Do not use prod near flammable liquids or dust.
- Place the safety clip under the trigger button whenever prod is not in use.
- Prod is intended for use on ground levels only.

Purpose

The SyrVet™ Livestock Prod is an electronic device designed to help move stubborn livestock.

Function

The SyrVet™ Livestock Prod contains a battery-powered electronic circuit which produces a high voltage pulse. The pulse is released by touching the animal with the two electrodes.

Animal Protection

Before using, refer to and follow the *Animal Protection and Slaughtering Policies and Guidelines* that pertain to your location.

Installing/Changing the Batteries

Unscrew the locking ring to separate the battery compartment from the motor/handle assembly. Remove the compartment door by turning it counter-clockwise. Install six new size “C” batteries as indicated by the diagrams on the inside of the battery compartment. Replace compartment door and turn it clockwise to lock into place. Replace locking ring and twist to reassemble motor/handle to battery compartment.

Testing the Prod

- 1) Install batteries as instructed above and as shown on the inside of the battery compartment.
- 2) Screw shaft onto the motor/handle assembly. Do not over-tighten.
- 3) Make sure the prod is put together correctly and all connections are secure.
- 4) Remove the safety clip from the trigger button.
- 5) Press the button and listen for the start of the motor (a buzzing sound).
- 6) Touch both prongs to a metal surface. Sparks should occur on contact.

Operating Instructions (Not to be used by children)

To operate the prod:

- 1) Install batteries as instructed above and as shown on the inside of the battery compartment.
- 2) Screw shaft onto the motor/handle assembly. Do not over-tighten.
- 3) Make sure the prod is put together correctly and all connections are secure.
- 4) Remove the safety clip from under the trigger button.
- 5) Touch the hind leg muscle of the animal with the two electrodes at the end of the prod.
- 6) Press and release the button. A light touch to the animal is the most effective. Only use the prod on an animal as necessary and according to the animal protection policies and guidelines of your area.

Maintenance & Storage: General Upkeep/Storage of the SyrVet™ Livestock Prod. Follow these maintenance guidelines after each use:

- 1) Completely CLEAN AND DRY the prod by detaching the shaft from handle and gently wiping the electrical connections on both the shaft and the handle with a dry cloth. Wipe the handle with a dry cloth, making sure to clean around the trigger button.
- 2) Examine the two electrode prongs at the shaft tip for any excessive wear or damage. Replace old shaft with a new one as necessary or if any signs of damage is visible.
- 3) When storing, separate shaft from handle to allow for complete drying in areas that can't be accessed. Only reassemble prod prior to use and only when it is thoroughly dry.
- 4) When not in use for long periods of time, remove batteries from the handle's battery compartment.
- 5) Store prod away from and out of reach of children.

Periodic Maintenance

Perform the following steps weekly or whenever prod has been exposed to moisture:

- 6) Disassemble prod. Clean and thoroughly dry all parts.
- 7) Examine and clean all metal contacts in battery compartment and at the shaft and motor/handle assembly connection for corrosion.
- 8) Test the battery power levels. Replace batteries if power level is weak or if batteries are corroded. SyrVet™ Livestock Prod Batteries are recommended. Replace all depleted batteries at the same time.
- 9) Examine all parts for damage. Replace parts as needed.

Replacement Parts

Replacement shafts, handles, batteries, motors, locking rings, and battery compartments for the SyrVet™ Livestock Prod are available through your local retailer or by calling SyrVet, Inc. at (800) 727-5203.

Troubleshooting


Follow these steps if the SyrVet™ Livestock Prod does not properly operate:

- 1) Check that the batteries are inserted properly and contain a charge. Do not use if batteries are corroded. Replace all depleted batteries at the same time. SyrVet™ Livestock Prod Alkaline Batteries are recommended for best performance.
- 2) Clean and examine all electrical contact points to make sure they are free of corrosion or damage. Replace worn, damaged, or old shafts, motors, or handles as needed.
- 3) Check that all parts are assembled properly and completely dry.

Visit www.syrvet.com for a more detailed troubleshooting guide.

Printed in China

SyrVet™ Livestock Prod

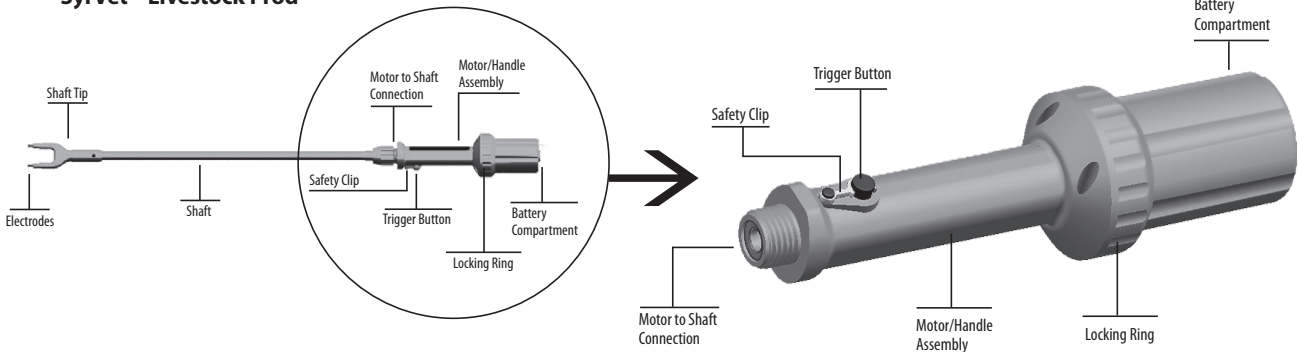SyrVet Livestock Prod Limited Warranty*

See reverse side for mail-in warranty form

SyrVet Livestock Prod Handle Unit—1 Year Limited Warranty* from date of purchase.

SyrVet Livestock Prod Shaft Only—Electrical continuity material and workmanship only. (*No warranty after sale to end user.*)

If the product fails due to a defect in the materials or workmanship within the warranty period, the product must be returned to SyrVet, Inc., postage paid for inspection and review. If the product is found to be defective, a replacement of the identical product will be returned.

**Limited Warranty is for prod electronics only. Shaft and/or handle breakage, any physical damage, or damages from normal wear and use is not covered. Warranty is non-transferable. To qualify for warranty service, the warranty registration information form along with a copy of the sales receipt must be mailed within 30 days of purchase to: SyrVet, Inc., Warranty Department, P.O. Box 490, Waukee, IA 50265*

ADDITIONALLY: To the maximum extent permitted by applicable law, the Company and its suppliers disclaim all other warranties, expressed or implied, including but not limited to implied warranties of merchantability and fitness for a particular purpose with regard to this product and any related written materials. Also, in no event shall the Company and its suppliers be liable for any damages whatsoever including without limitation special, incidental, consequential, or indirect damages for personal injury, loss of business profits, or any other pecuniary losses arising out of the use or inability to use this product. The Company and its supplier's entire liability shall be limited to the amount actually paid by the consumer for this product. Because some State jurisdictions do not allow the exclusion or limitation of liability for consequential or incidental damages, the above limitations may not apply to you.

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SYRVET LIMITED WARRANTY REGISTRATION FORM

Thank you for purchasing your SyrVet Livestock Prod!

To qualify for the Limited Warranty as described on the reverse side, please complete the following information and mail to the address below.

NAME: _____ DATE OF PURCHASE: _____

STREET ADDRESS/CITY: _____ STATE: _____ ZIP CODE: _____

E-MAIL ADDRESS: _____ PHONE #: _____

PURCHASED FROM: _____ CITY/STATE: _____

*Please attach a copy of the receipt from the purchase of your SyrVet Livestock Prod along with this warranty registration form within 30 days from date of purchase to: **SyrVet, Inc., Warranty Dept., P.O. Box 490, Waukee, IA 50265***

Make a copy of this form (front and back) for your records.